

QA/QC Compliance Assessment to assist with Quality Review

Work Order : **EM2115784** Page : 1 of 5

Client : Dept for Environment & Water : Laboratory : Environmental Division Melbourne

 Contact
 : Mr FRANK MANGERUCA
 Telephone
 : +61881625130

 Project
 : HCHB
 Date Samples Received
 : 11-Aug-2021

 Site
 : --- Issue Date
 : 16-Aug-2021

Sampler : JC No. of samples received : 2
Order number : ---- No. of samples analysed : 2

This report is automatically generated by the ALS LIMS through interpretation of the ALS Quality Control Report and several Quality Assurance parameters measured by ALS. This automated reporting highlights any non-conformances, facilitates faster and more accurate data validation and is designed to assist internal expert and external Auditor review. Many components of this report contribute to the overall DQO assessment and reporting for guideline compliance.

Brief method summaries and references are also provided to assist in traceability.

Summary of Outliers

Outliers: Quality Control Samples

This report highlights outliers flagged in the Quality Control (QC) Report.

- NO Method Blank value outliers occur.
- NO Duplicate outliers occur.
- NO Laboratory Control outliers occur.
- NO Matrix Spike outliers occur.
- For all regular sample matrices, NO surrogate recovery outliers occur.

Outliers: Analysis Holding Time Compliance

• Analysis Holding Time Outliers exist - please see following pages for full details.

Outliers : Frequency of Quality Control Samples

• NO Quality Control Sample Frequency Outliers exist.


Page 2 of 5 EM2115784 Work Order

Dept for Environment & Water Client

Project



Matrix: WATER

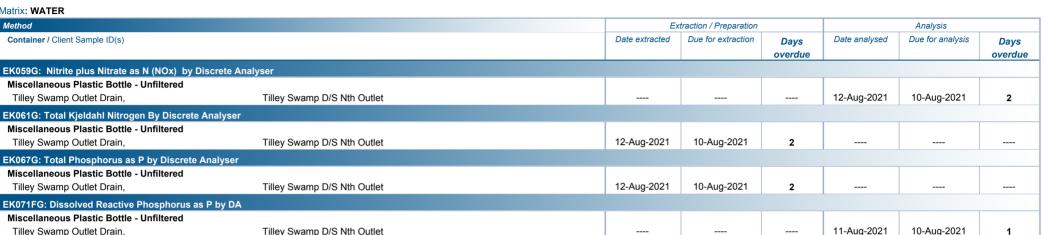

Analysis Holding Time Compliance

If samples are identified below as having been analysed or extracted outside of recommended holding times, this should be taken into consideration when interpreting results,

This report summarizes extraction / preparation and analysis times and compares each with ALS recommended holding times (referencing USEPA SW 846, APHA, AS and NEPM) based on the sample container provided. Dates reported represent first date of extraction or analysis and preclude subsequent dilutions and reruns. A listing of breaches (if any) is provided herein.

Holding time for leachate methods (e.g. TCLP) vary according to the analytes reported. Assessment compares the leach date with the shortest analyte holding time for the equivalent soil method. These are: organics 14 days, mercury 28 days & other metals 180 days. A recorded breach does not quarantee a breach for all non-volatile parameters.

Holding times for VOC in soils vary according to analytes of interest. Vinyl Chloride and Styrene holding time is 7 days; others 14 days. A recorded breach does not guarantee a breach for all VOC analytes and should be verified in case the reported breach is a false positive or Vinyl Chloride and Styrene are not key analytes of interest/concern.

Matrix: WATER

Evaluation: x = Holding time breach: $\sqrt{}$ = Within holding time

IVIAUIA. WATER					Lvaluation	. ~ - Holding time	breach, • - with	ir noluling tilli
Method			Extraction / Preparation			Analysis		
Container / Client Sample ID(s)		Date extracted	Due for extraction	Evaluation	Date analysed	Due for analysis	Evaluation	
EA015: Total Dissolved Solids dried at 180 ± 5 °C								
Miscellaneous Plastic Bottle - Unfiltered (EA015H) Tilley Swamp Outlet Drain,	Tilley Swamp D/S Nth Outlet	09-Aug-2021				13-Aug-2021	16-Aug-2021	✓
EG052G: Silica by Discrete Analyser								
Miscellaneous Plastic Bottle - Unfiltered (EG052G) Tilley Swamp Outlet Drain,	09-Aug-2021				13-Aug-2021	06-Sep-2021	✓	
EK055G-SW: Ammonia as N by Discrete Analyser in S	aline Water							
Miscellaneous Plastic Bottle - Unfiltered (EK055G-SW) Tilley Swamp Outlet Drain,	Tilley Swamp D/S Nth Outlet	09-Aug-2021				16-Aug-2021	06-Sep-2021	✓
EK057G: Nitrite as N by Discrete Analyser								
Miscellaneous Plastic Bottle - Unfiltered (EK057G) Tilley Swamp Outlet Drain,	Tilley Swamp D/S Nth Outlet	09-Aug-2021				11-Aug-2021	11-Aug-2021	✓

Page : 3 of 5
Work Order : EM2115784

Client : Dept for Environment & Water

Project : HCHB



Matrix: WATER					Evaluation	: × = Holding time	breach ; ✓ = Withi	n holding time
Method			Extraction / Preparation			Analysis		
Container / Client Sample ID(s)			Date extracted	Due for extraction	Evaluation	Date analysed	Due for analysis	Evaluation
EK059G: Nitrite plus Nitrate as N (NOx) by Discrete Ana	llyser							
Miscellaneous Plastic Bottle - Unfiltered (EK059G) Tilley Swamp Outlet Drain,	Tilley Swamp D/S Nth Outlet	09-Aug-2021				12-Aug-2021	10-Aug-2021	*
EK061G: Total Kjeldahl Nitrogen By Discrete Analyser								
Miscellaneous Plastic Bottle - Unfiltered (EK061G) Tilley Swamp Outlet Drain,	Tilley Swamp D/S Nth Outlet	09-Aug-2021	12-Aug-2021	10-Aug-2021	<u>se</u>	13-Aug-2021	09-Sep-2021	1
EK067G: Total Phosphorus as P by Discrete Analyser								
Miscellaneous Plastic Bottle - Unfiltered (EK067G) Tilley Swamp Outlet Drain,	Tilley Swamp D/S Nth Outlet	09-Aug-2021	12-Aug-2021	10-Aug-2021	<u>Je</u>	13-Aug-2021	09-Sep-2021	✓
EK071FG: Dissolved Reactive Phosphorus as P by DA								
Miscellaneous Plastic Bottle - Unfiltered (EK071FG) Tilley Swamp Outlet Drain,	Tilley Swamp D/S Nth Outlet	09-Aug-2021				11-Aug-2021	10-Aug-2021	JC

Page : 4 of 5 Work Order : EM2115784

Client : Dept for Environment & Water

Project : HCHI



Quality Control Parameter Frequency Compliance

The following report summarises the frequency of laboratory QC samples analysed within the analytical lot(s) in which the submitted sample(s) was(were) processed. Actual rate should be greater than or equal to the expected rate. A listing of breaches is provided in the Summary of Outliers.

	<u> </u>						not within specification; ✓ = Quality Control frequency within specific
Quality Control Sample Type	14 (/ /		ount		Rate (%)	Fratratia	Quality Control Specification
Analytical Methods	Method	QC	Reaular	Actual	Expected	Evaluation	
aboratory Duplicates (DUP)							
Ammonia as N (Saline Water)	EK055G-SW	1	8	12.50	10.00	✓	NEPM 2013 B3 & ALS QC Standard
Dissolved Reactive Phosphorus as P by DA	EK071FG	1	2	50.00	10.00	✓	NEPM 2013 B3 & ALS QC Standard
Nitrite and Nitrate as N (NOx) by Discrete Analyser	EK059G	1	8	12.50	10.00	✓	NEPM 2013 B3 & ALS QC Standard
litrite as N by Discrete Analyser	EK057G	2	18	11.11	10.00	✓	NEPM 2013 B3 & ALS QC Standard
filica (Reactive) by Discrete Analyser	EG052G	1	3	33.33	10.00	✓	NEPM 2013 B3 & ALS QC Standard
otal Dissolved Solids (High Level)	EA015H	2	14	14.29	10.00	✓	NEPM 2013 B3 & ALS QC Standard
otal Kjeldahl Nitrogen as N By Discrete Analyser	EK061G	2	18	11.11	10.00	✓	NEPM 2013 B3 & ALS QC Standard
otal Phosphorus as P By Discrete Analyser	EK067G	2	15	13.33	10.00	✓	NEPM 2013 B3 & ALS QC Standard
aboratory Control Samples (LCS)							
Ammonia as N (Saline Water)	EK055G-SW	1	8	12.50	5.00	✓	NEPM 2013 B3 & ALS QC Standard
Dissolved Reactive Phosphorus as P by DA	EK071FG	1	2	50.00	5.00	✓	NEPM 2013 B3 & ALS QC Standard
litrite and Nitrate as N (NOx) by Discrete Analyser	EK059G	1	8	12.50	5.00	✓	NEPM 2013 B3 & ALS QC Standard
itrite as N by Discrete Analyser	EK057G	1	18	5.56	5.00	√	NEPM 2013 B3 & ALS QC Standard
ilica (Reactive) by Discrete Analyser	EG052G	1	3	33.33	5.00	✓	NEPM 2013 B3 & ALS QC Standard
otal Dissolved Solids (High Level)	EA015H	2	14	14.29	10.00	√	NEPM 2013 B3 & ALS QC Standard
otal Kjeldahl Nitrogen as N By Discrete Analyser	EK061G	1	18	5.56	5.00	√	NEPM 2013 B3 & ALS QC Standard
otal Phosphorus as P By Discrete Analyser	EK067G	1	15	6.67	5.00	✓	NEPM 2013 B3 & ALS QC Standard
lethod Blanks (MB)							
mmonia as N (Saline Water)	EK055G-SW	1	8	12.50	5.00	1	NEPM 2013 B3 & ALS QC Standard
Dissolved Reactive Phosphorus as P by DA	EK071FG	1	2	50.00	5.00	√	NEPM 2013 B3 & ALS QC Standard
litrite and Nitrate as N (NOx) by Discrete Analyser	EK059G	1	8	12.50	5.00	√	NEPM 2013 B3 & ALS QC Standard
litrite as N by Discrete Analyser	EK057G	1	18	5.56	5.00	√	NEPM 2013 B3 & ALS QC Standard
ilica (Reactive) by Discrete Analyser	EG052G	1	3	33.33	5.00	√	NEPM 2013 B3 & ALS QC Standard
otal Dissolved Solids (High Level)	EA015H	1	14	7.14	5.00	√	NEPM 2013 B3 & ALS QC Standard
otal Kjeldahl Nitrogen as N By Discrete Analyser	EK061G	1	18	5.56	5.00	√	NEPM 2013 B3 & ALS QC Standard
otal Phosphorus as P By Discrete Analyser	EK067G	1	15	6.67	5.00	√	NEPM 2013 B3 & ALS QC Standard
Matrix Spikes (MS)							
mmonia as N (Saline Water)	EK055G-SW	1	8	12.50	5.00	✓	NEPM 2013 B3 & ALS QC Standard
Dissolved Reactive Phosphorus as P by DA	EK071FG	1	2	50.00	5.00		NEPM 2013 B3 & ALS QC Standard
litrite and Nitrate as N (NOx) by Discrete Analyser	EK059G	1	8	12.50	5.00	<u> </u>	NEPM 2013 B3 & ALS QC Standard
litrite as N by Discrete Analyser	EK057G	1	18	5.56	5.00		NEPM 2013 B3 & ALS QC Standard
ilica (Reactive) by Discrete Analyser	EG052G	1	3	33.33	5.00		NEPM 2013 B3 & ALS QC Standard
otal Kjeldahl Nitrogen as N By Discrete Analyser	EK061G	1	18	5.56	5.00	<u> </u>	NEPM 2013 B3 & ALS QC Standard
otal Phosphorus as P By Discrete Analyser	EK067G	1	15	6.67	5.00	<u> </u>	NEPM 2013 B3 & ALS QC Standard

Page : 5 of 5 Work Order : EM2115784

Client : Dept for Environment & Water

Project : HCHI



Brief Method Summaries

The analytical procedures used by the Environmental Division have been developed from established internationally recognized procedures such as those published by the US EPA, APHA, AS and NEPM. In house developed procedures are employed in the absence of documented standards or by client request. The following report provides brief descriptions of the analytical procedures employed for results reported in the Certificate of Analysis. Sources from which ALS methods have been developed are provided within the Method Descriptions.

Analytical Methods	Method	Matrix	Method Descriptions
Total Dissolved Solids (High Level)	EA015H	WATER	In house: Referenced to APHA 2540C. A gravimetric procedure that determines the amount of `filterable` residue in an aqueous sample. A well-mixed sample is filtered through a glass fibre filter (1.2um). The filtrate is evaporated to dryness and dried to constant weight at 180+/-5C. This method is compliant with NEPM Schedule B(3)
Silica (Reactive) by Discrete Analyser	EG052G	WATER	In house: Referenced to APHA 4500-SiO2 D: Under Acdic conditions reactive silicon combines with ammonium molybdate to form a yellow molybdosilicic acid complex. This is reduced by 1-amino-2-naphthol-4-sulfonic acid to a silicomolybdenum blue complex which is measured by discrete analyser at 670 nm. This method is compliant with NEPM Schedule B(3).
Ammonia as N (Saline Water)	EK055G-SW	WATER	In house: Referenced to APHA 4500-NH3 G Ammonia is determined by direct colorimetry by Discrete Analyser. This method is compliant with NEPM Schedule B(3)
Nitrite as N by Discrete Analyser	EK057G	WATER	In house: Referenced to APHA 4500-NO2- B. Nitrite is determined by direct colourimetry by Discrete Analyser. This method is compliant with NEPM Schedule B(3)
Nitrate as N by Discrete Analyser	EK058G	WATER	In house: Referenced to APHA 4500-NO3- F. Nitrate is reduced to nitrite by way of a chemical reduction followed by quantification by Discrete Analyser. Nitrite is determined seperately by direct colourimetry and result for Nitrate calculated as the difference between the two results. This method is compliant with NEPM Schedule B(3)
Nitrite and Nitrate as N (NOx) by Discrete Analyser	EK059G	WATER	In house: Referenced to APHA 4500-NO3- F. Combined oxidised Nitrogen (NO2+NO3) is determined by Chemical Reduction and direct colourimetry by Discrete Analyser. This method is compliant with NEPM Schedule B(3)
Total Kjeldahl Nitrogen as N By Discrete Analyser	EK061G	WATER	In house: Referenced to APHA 4500-Norg D (In house). An aliquot of sample is digested using a high temperature Kjeldahl digestion to convert nitrogenous compounds to ammonia. Ammonia is determined colorimetrically by discrete analyser. This method is compliant with NEPM Schedule B(3)
Total Nitrogen as N (TKN + Nox) By Discrete Analyser	EK062G	WATER	In house: Referenced to APHA 4500-Norg / 4500-NO3 This method is compliant with NEPM Schedule B(3)
Total Phosphorus as P By Discrete Analyser	EK067G	WATER	In house: Referenced to APHA 4500-P H, Jirka et al, Zhang et al. This procedure involves sulphuric acid digestion of a sample aliquot to break phosphorus down to orthophosphate. The orthophosphate reacts with ammonium molybdate and antimony potassium tartrate to form a complex which is then reduced and its concentration measured at 880nm using discrete analyser. This method is compliant with NEPM Schedule B(3)
Dissolved Reactive Phosphorus as P by DA	EK071FG	WATER	In house: Referenced to APHA 4500-P F Water samples are filtered through a 0.45um filter prior to analysis. Ammonium molybdate and potassium antimonyl tartrate reacts in acid medium with othophosphate to form a heteropoly acid -phosphomolybdic acid - which is reduced to intensely coloured molybdenum blue by ascorbic acid. Quantification is achieved by Discrete Analyser. This method is compliant with NEPM Schedule B(3)
Preparation Methods	Method	Matrix	Method Descriptions
TKN/TP Digestion	EK061/EK067	WATER	In house: Referenced to APHA 4500 Norg - D; APHA 4500 P - H. This method is compliant with NEPM Schedule B(3)